# SONIQ

## SONIQ Lite

Dentium  
For Dentists By Dentists

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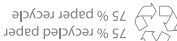
### Component

English

\* SONIQ Lite package is not include 'Regular head(3PCS)' and 'Travel Case'.

Toothbrush product

Travel case

Wireless charging cradle

Regular head(3PCS)

User manual

1

### Product structure

English

Toothbrush head

Power/Action button

Mode indicator LED


Charging indicator LED

Wireless charging cradle


2

### Product function

English
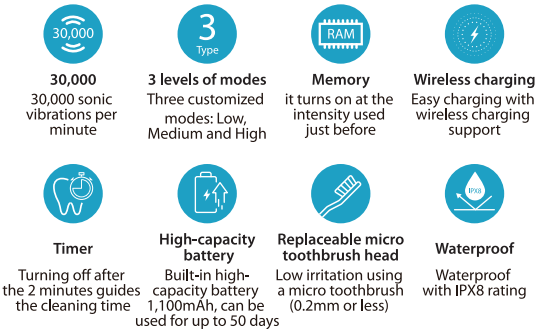
#### Toothbrush head replacement

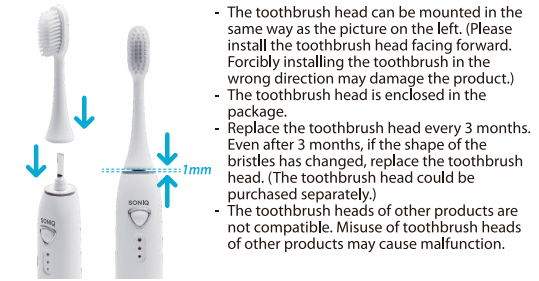- The toothbrush head can be mounted in the  
same way as the picture on the left. (Please  
install the toothbrush head facing forward.  
Forcibly installing the toothbrush in the  
wrong direction may damage the product.)
- The toothbrush head is enclosed in the  
package.
- Replace the toothbrush head every 3 months.  
Even after 3 months, if the shape of the  
bristles has changed, replace the toothbrush  
head. (The toothbrush head could be  
purchased separately.)
- The toothbrush heads of other products are  
not compatible. Misuse of toothbrush heads  
of other products may cause malfunction.

3


English

- After combining the toothbrush head, there is a gap of about 1mm between the body and the toothbrush head. This gap is a space for the normal work of the brush head. (If you press excessively to close the gap, it may cause the product failure.)

#### Usage Method

- 1. Prepare to brush your teeth**  
Take a small amount of toothpaste on the brush head. Vibration increases gradually for 3 seconds and this prevents the toothpaste on the brush splashed. Put it into the mouth within 3 seconds.
- 2. Power ON, mode control**  
Briefly press the power button to turn on. While the power is on, briefly press the power button again to change the mode and start brushing on the mode you prefer.
- 3. Use, power OFF**  
We recommend brushing your teeth for about 2 minutes. The power turns off automatically after 2 minutes. If you want to turn off the power while brushing your teeth, press and hold the power button for more than 1 second.
- 4. Wash and dry**  
After use, rinse the toothbrush head under running water and let it dry naturally. After use, store in a dry and cool place.

#### Smart 2 minutes Timer: 30 Seconds Interval

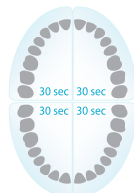- This product is equipped with a timer to  
guide the brushing time.
- Brushing your teeth for 2 minutes and turns  
off automatically.
- Brush your teeth in each zone of the oral  
cavity and change the zone when it vibrates  
every 30 seconds.

4